



## COMPARE MODELS

	3100A	3000A	3200A
Sample Capacity	15 samples	121 samples (1 removable rack)	209 samples (2 removable racks)
User Interface	Keypad	Touch Screen	Touch Screen
Easy fixing/ Repositioning system	✓	✓	✓
AI functionalities	✓	✓	✓
Virtual Touch Screen/ Screen Mirroring	✓	✓	✓
Syringe Housing Illumination	-	✓	✓
SyringeID <sup>2</sup>	-	-	✓
BCR (Bar Code Reader) <sup>1</sup>	-	-	Optional

## TECHNICAL SPECIFICATIONS

### General features

Syringe volume<sup>1</sup>: 0.5, 1, 1.2, 5, 10, 25, 50 and 100µl  
 Maintenance: preventive maintenance capabilities  
 Electrical control: LAN and TTL; optional: RS232  
 Syringe area illumination: yes (programmable) for 3000A/3200A

### Tray capacity<sup>4</sup>

3000A: 121 vials, 2ml (1 removable rack)  
 3200A: 209 vials, 2ml (2 removable racks)  
 3100A: 15 vials, 2ml  
 A thermally conditioned rack is available for 3000A as accessory

### Filling

Sample volume: as low as step of 0.1µl  
 Air volume: as low as step of 0.1µl  
 Filling speed: 1-100µl/sec  
 Viscosity delay: 0-15s  
 Bubble elimination: up to 15 pull up strokes

### Injection

Injection speed: 1-100µl/sec  
 Injection depth: programmable  
 Pre and post inj delay: 0-99s

### Washing

Type: pre-injection, sample, post-injection  
 Solvent capacity: 6x10ml vials  
 Mode: single or double wash

### Internal standard technique

IS volume: as low as step of 0.1µl  
 Air gap volume: as low as step of 0.1µl  
 Mode: 1 or 2 air gaps

### Physical features

Dimensions (WxHxD)<sup>5</sup>: 280x570x320mm  
 Weight (3000A): 8.0kg  
 Weight (3200A): 9.2kg  
 Weight (3100A): 6.4kg  
 Power supply: 100-240±10%Vac; 50-60Hz; 55W

### Software

HTA Monitor: included for free  
 Virtual Touch Screen and AI functionalities  
 60-days free trial  
 Full autosampler PC programming

### HTA Monitor PC requirements

Software:
 

- Microsoft Windows 7, Windows 8.1, Windows 10, Windows 11 PC Edition only (excluding mobile devices and appliances)
- PC is expected to run Windows OS with the latest update installed (unless differently specified)
- Additional required software: Microsoft .NET Framework 4.5.2
- Administrator account for software installation only

 Hardware:
 

- RAM: 2GB
- Disk space (for installation): 6GB
- LAN port
- 1024x768 Minimum
- Some functionalities require Internet access

<sup>1</sup> Not compatible with Serial Port (RS232)

<sup>2</sup> Patented technology

<sup>4</sup> OEM and Custom versions with extended syringe range and custom trays are available

<sup>5</sup> Tray in closed position for 3000A and 3200A

<sup>6</sup> Subscription required. Initial 3-year subscription is provided with every instrument

HTA Monitor and Internet connection are required for some AI functionalities.  
 It is not required to have HTA Monitor and CDS on the same PC.  
 PC is required for setup, service and programming 3100A.  
 All trademarks are property of their respective owners.



At HTA, we design and manufacture robotics solutions for the scientific community. With decades of expertise, we deliver reliable technologies tailored to modern laboratories. Our products for analytical chemistry include autosamplers and sample prep solutions, covering, among others, GC, LC, and ICP apps. We offer universal autosamplers compatible with analyzers from multiple brands, ensuring flexible lab workflows. HTA products are made in Italy under UNI EN ISO 9001:2015 and 13485:2021 certified quality systems.

### HTA s.r.l.

via del Mella, 21 - 25131 Brescia - ITALY  
 T: +39 030 3582920  
[www.hta-it.com](http://www.hta-it.com) | [enquiry@hta-it.com](mailto:enquiry@hta-it.com)



Distributed by:



AI POWERED



## 3000A series LIQUID GC AUTOSAMPLERS


Offering reliability and superior performance for all GC-GC/MS applications, meeting the requirements of both today and tomorrow.

## 3200A

Due to its 209 2mL sample vials, the **3200A** has set a new standard in the market. Samples are organised in two racks that can be easily removed for sample loading or preparation, or to be stored elsewhere.

An **integrated Bar Code Reader**<sup>1</sup> can be available with every **3200A** unit, offering full and detailed sample tracking.

The **3200A** features **SyringeID**<sup>2</sup>, a patented technology based on RFID tags. The SyringeID can **identify syringes** in a univocal way; thereby preventing errors when mounting a syringe, preventing syringe volume mismatching and keeps track of the syringe consumption (preventive maintenance counters). The SyringeID system can provide you with a level of confidence never previously achieved.

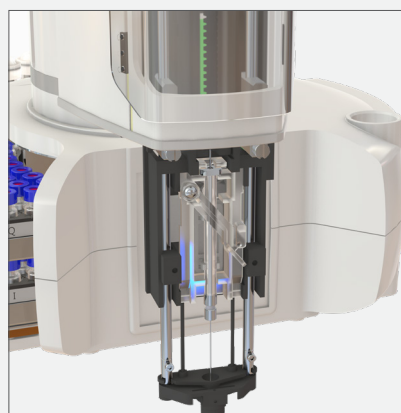

## KEY FEATURES:

- 60+ GC models supported
- Space-saving design
- Move freely between GCs - 1 turret serves 2 injectors
- Extraordinary performance: a perfect fit for your lab
- Empowered by AI



The **3000A series** has been engineered for **maximum performance and reliability**. It takes advantage of today's latest technology to deliver even greater reliability, performance and flexibility.

## SPACE-SAVING DESIGN

**3000A series** autosamplers are the **most compact** in the market, with a near-to-zero requirement for bench space, while offering a **top-class sample capacity**.

The **rotating tower leaves the injector port free for manual injection or maintenance**. In addition to avoiding any complicated sample vial transfer procedures, this approach keeps the **samples away from any heated source**. The sample rack is mounted away from the GC oven to prevent exposure to high temperatures, which could cause degradation or condensation in the sample vial.

The **system is fully self-contained and can be easily installed on all the GC and GC/MS systems** available in the market, because of its flexibility and modular configuration. It can be configured to **serve up to two injectors** in most supported GCs.

## EASY TO OPERATE

The self-aligning "plug and play" autosamplers of the **3000A series** mount in seconds and without the need for tools. It can be easily moved between GCs when workloads change, due to its easy positioning, repositioning and easy removal.

**Easy to use:** Just load the samples and run the analysis with no extra

downtime. **Run-all-samples mode** makes the operations of sample processing simple and transparent; the sample rack can be removed for sample loading or preparation or stored elsewhere.

The **high-quality touchscreen** provides easier system accessibility and usability: it eliminates drill down, simplifying instrument control for both novices and experienced users. Besides the touch screen, the **3000A series** can be also controlled by a PC using the HTA Autosampler Manager (please see the dedicated brochure), available in standard or CFR 21 Part 11 compliant version.

**Syringe illumination** simplifies method validation by allowing to keep the sample under control while a **quick-fix syringe bracket** allows easy replacement of the syringe. Prevention of errors through **vial detection**: **3000A series** autosamplers can detect the presence or absence of vials for sample, solvent and waste, to notify the analyst or skip to the next sample if programmed so.

## A PERFECT FIT FOR YOUR LAB

With HTA, you have the flexibility to choose the automatic sampler that best suits the needs of your lab. The **3000A series** is available in the standard version with 121 positions for 2mL vials. The top version **3200A** has 209 positions, while the entry-level version **3100A** has a 15-position fixed rack.

The **HTA GC autosamplers are the fastest liquid autosamplers available** on the market: injection is performed in less than 100ms. Fast-injection technology ensures the best possible peak shape, while maximising the accuracy of your results.

The **great solvent capacity** means longer unattended operations: six solvent vials of 10mL each, giving a total capacity of 60mL.

The ability to customise cleaning solvent types and sequences includes the **double wash step capability**: as pre- and post-wash solvent in addition to A, B... F solvents you can also choose a combination of A+B, A+C... F+E solvents for superior analytical performance (carry-over adverse!).

The **3000A series autosamplers** manage the **most sophisticated sampling techniques**: internal standard technique (also known as sandwich injection), multiphase, ambient headspace, priority injection, nanoliter and large volume injection. Parameters are easily programmable to optimise both the most convenient sampling methods for both extremely volatile or viscous samples and the best injection technique. **Variable needle depths** allow samples to be taken from anywhere within the vial, thus performing extractions directly from the vial.

## EMPOWER YOUR MASS SPECTROMETER

Make the best choice to address evolving analytical needs: specific functionalities have been engineered to elevate MS potentialities and next-generation analysers.

Mass detectors and next-gen analysers are more sensitive than a normal GC to some phenomena: that's the reason the 3000A series **exact and highly reproducible position of the syringe needle in the injector port** makes the difference. Furthermore, **sampling and injection methods reducing stress to septa** have been implemented: this leads to minimising contamination of liner and analyser.


Special functionalities to support prep-run have been implemented to **reduce carrier gas consumption** on the analyser resulting in lower operations cost, increased analyser parts lifetime and a more environmentally-conscious behaviour.

## 3100A

Analysing fewer samples? Get precision and performance with a system that is perfectly sized for your needs. The **3100A** offers automation for processing up to 15 unattended samples.

Is limited interaction preferred? A **single button start-stop approach** is the best solution: quick sets of autosampler parameters can be done by PC, while all the routine operations such as analysis start, sample loading and extra wash can be managed by its dedicated keypad.

The **new led bar and PC quick link makes everything smooth**: color-coded status indicator provides at-a-glance autosampler status information that is easy to see; virtual screen on PC shows up automatically when an informative message is available; alternatively, it can be activated by double thumb pressure on keypad.

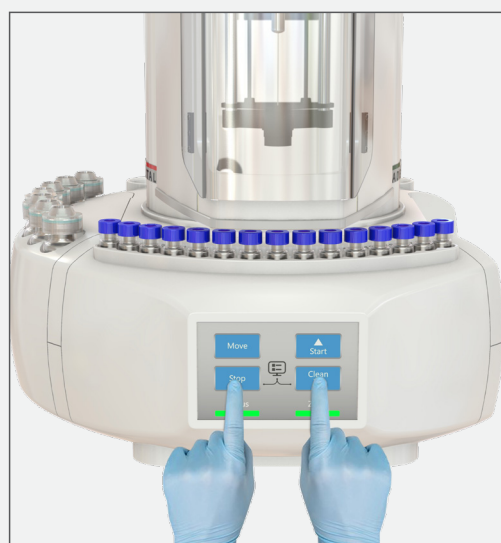

# ARTIFICIAL INTELLIGENCE

HTA Monitor – PC utility – is the engine at the foundation of our Artificial Intelligence (AI) capabilities<sup>3</sup>. Learn below how AI can boost your lab productivity!



## GC INLET SEPTUM TEST

Concerned about having an overtightened GC inlet septum nut after **septum replacement**? Don't worry: the **3000A series** provides the ability to **check whether the septum is overtightened** or correctly compressed. The excess of septum compression can result in fact in septum corning (reduced septum lifetime, fragments of septum introduced into the injection port) or – in extremely cases – to needle damage (bent needle or shorten needle life).



## AUTOMATED CONSUMABLES TRACKING

Automated consumables tracking with notifications **minimises unexpected downtime** and waste due to unnecessary replacement. Consumables consumption tracking goes far above preventive maintenance counters: **expiration dates and performance tests** are available. On each start-up, a **syringe plunger diagnostic test** is performed to check whether syringe plunger maintenance or syringe replacement is needed.



## SUPPORT FLEXIBLE WORKING STYLES

**Virtual screen (3100A)** and **screen mirroring (3000A and 3200A)** enable you to control the autosampler from your PC without the need to stand in front of the autosampler. A replica of your autosampler touch screen is made available to you so that you can perform every task from the same familiar user interface. Our products have an **automated email alert system**<sup>4</sup>: in case of an interruption, an email is sent instantly, enabling quick action to minimize downtime and avoid delays.



## EASY CONNECT WITH SERVICE

You can **contact tech support by scanning a QR code** and passing all relevant information about your instrument, configuration and issue. In most cases, before you even ask a question, we will have given you the answer!



## EMF+: ADVANCED EARLY MAINTENANCE FEEDBACK

While traditional preventive maintenance ensures robustness by "over-maintaining" your equipment, EMF+ shifts the paradigm toward true intelligence. Our **AI engine** goes beyond simple counters, performing **automated self-diagnostic tests** during idle times or upon request. By analyzing the actual instrument status, EMF+ predicts maintenance needs with surgical precision. The result? The same rock-solid reliability of preventive maintenance, but at a fraction of the cost and with significantly extended uptime.



## INSTRUMENT PARAMETER OPTIMIZATION AND LIVE INSTRUCTIONS

For key tasks, our **live instructions provide video guides** to refresh experienced operators and train new staff effectively. Additionally, they enhance the user experience by **suggesting programming and setting adjustments** for a smoother, more efficient analysis flow.

The **AI engine is continuously improving as it regularly receives software updates via the internet**. These updates add new functionalities and improve existing ones: enable the automatic update function to always stay up to date. AI engine updates do not affect the operations of the autosamplers and this is why they are safe and well-accepted even in highly regulated contexts.